


BTS2048 Series

<https://www.gigahertz-optik.com/en-us/product/bts2048-series/>

Product tags: UV , VIS , NIR , Weatherproof

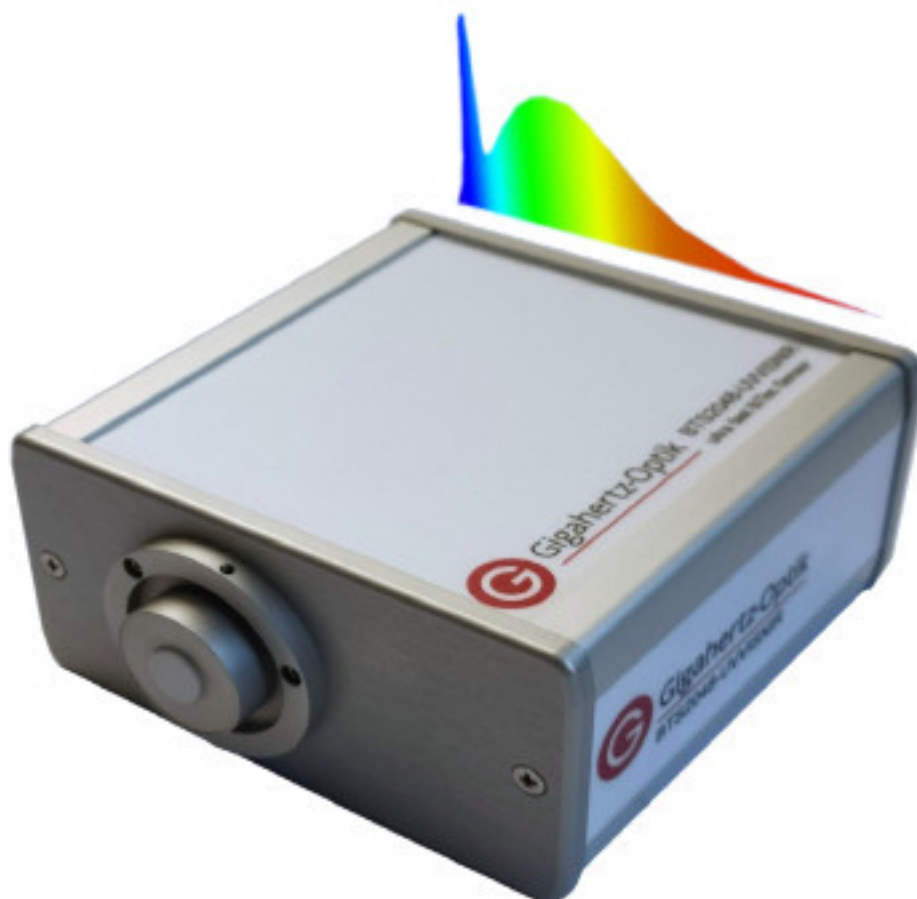

Description

Meets the Requirements of a Modern, Versatile Array Spectroradiometer from UV to IR

There are many factors to be considered when assessing the suitability of array spectrometer systems for the absolute measurement of optical radiation parameters.

- For instance, the measurement of lamps that have different power ratings is only possible using array detectors that have a wide dynamic range. Precise, absolute measurements require the entire dynamic range of the spectroradiometer to be completely linearized and also require an accurate, traceable calibration.
- If the electronically controlled dynamic range (set by the integration time) is not sufficient, additional attenuating filters are needed. The filter selector must be mechanically robust to ensure long-term stability of the measurement system. For time-critical applications such as LED binning in pulsed mode, the electronically controlled dynamic range must be large enough to avoid the need for a time-consuming filter change during the measurement. For absolute measurements, an automated dark signal adjustment of the CCD is most beneficial.
- Spectroradiometers that are used for the binning of LEDs must accommodate the precise synchronization of the measurement with the test LEDs operated in pulsed mode, thereby requiring suitable trigger interface and fast data readout. Flash measurements, i.e. measurements within a light pulse, require an electronic shutter for instantaneous (ns) zero setting of all pixels before a measurement is triggered.
- For accurate measurements in the UV range stray light reduction comparable to that of a double-monochromator is needed. Here the performance of the BTS2048 series is outstanding due to its combined optical filter and stray light matrix correction approach. This exceptional performance of the BTS2048 series is evidenced in several [scientific publications](#).
- For low intensity detection, long integration times are needed which require precise temperature stabilized sensors.
- For outdoor applications weatherproof housings are needed.
- For flicker measurements a fast CW integral sensor is needed. The [BTS technology](#) incorporates such an integral detector, which can be also be used as a reference sensor, thereby making the BTS2048 series superior to simple spectroradiometers.
- For low measurement uncertainties a high linearity (f_3), low f_2 error and low calibration uncertainties are needed. Our [ISO 17025 calibration laboratory](#) is the ideal facility for achieving this.
- The measurement of the luminous flux, luminous intensity, and luminous intensity distribution require additional accessory components e.g., integrating spheres, luminous intensity lenses, fiber coupled adapters and goniometers. Reproducible interfacing to these entrance optic accessories is essential. Direct mounting of the spectroradiometer onto the accessory equipment helps avoid influences of flexible light guide connections.
- Among all the above requirements, system performance in accordance with CIE 15, and TM-30-18, CIE S025 and LM-79-08 standards, etc. is needed. We are active partner in the DIN and CIE to name a few examples. This allows us to be up to date with latest [certification](#).
- For greatest versatility a compact size and highest mechanical robustness are essential.

The BTS2048 series combines this performance requirements throughout the product range. For almost every application a version of the BTS2048 is available. If not, a custom version is also an option.

Scientific Grade, Industrial Standard

The BTS2048 series has shown its performance in many scientific publications including inter-comparisons with double monochromators in terms of stray light reduction (see [scientific articles](#)). Furthermore, the well established BTS technology offers the possibility to have a reference sensor implemented, a fast diode for flicker measurements or simply a very stable integral detector. Additionally the integral diode can be spectral mismatch corrected by the spectral data and therefore is able to correct the spectral data in linearity/stability. See all advantages in this [technical article](#).

Due to this performance the BTS2048 series is applied world wide in LED measurement applications, as a UV reference laboratory device, for risk group determination of light sources, as a spectral irradiance meter in goniometers and many more.



Integrating Sphere: Direct mounting of the measurement device and accessory components



Fiber coupled Version for Gigahertz Optik fibers or FC



Compact size and orientation independent use for goniometric measurements

▢ See our detailed product series description (tab [Downloads](#))

Product Overview

[Here](#) you will find a product overview according to the applications.

In the following table the main versions configured for spectral irradiance are listed:

| Series | Spectral Range / nm |
|---|---------------------|
| BTS2048-UV Series | 200 - 430 |
| BTS2048-UV-2 Series | 200 - 550 |
| BTS2048-BS Series | 400 - 530 |
| BTS2048-UVVISNIR Series | 200 - 900 |
| BTS2048-VL-TEC Series | 280 - 1050 |
| BTS2048-IR Series | 950 - 2150 |

User Software and Software Development Kit

The standard S-BTS2048 user software offers a customizable user interface and is intuitive and easy to use. It has a large number of display and function modules which can be activated when configuring the BTS2048-series with the respective accessories from Gigahertz Optik GmbH. Furthermore customer specific software tools are possible if needed.

The S-SDK-BTS2048 developer software is recommended for integration of the BTS2048-series in the customer's own software.

Calibration


One essential quality feature of spectroradiometric devices is their precise and traceable calibration. The BTS2048-series is calibrated by Gigahertz-Optik's calibration laboratory that was accredited by DAkkS (D-K-15047-01-00) for the *spectral responsivity* and *spectral irradiance* according to ISO/IEC 17025. The calibration also includes the corresponding accessory components. Every device is delivered with its respective calibration certificate.



Weather Proof Housing for Outdoor Usage (UV Index, TOC, spectral irradiance, global UV)



SST-180x System Configuration 200 nm bis 3000 nm. Adapter plate for tripod available.

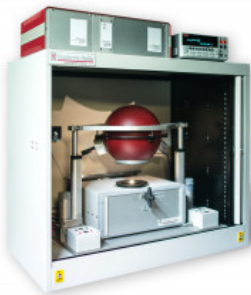


LED characterication Test Rig



BTS2048-VL-F-Z01: Probe header for small LED chips (waferprober)

THE BENEFITS

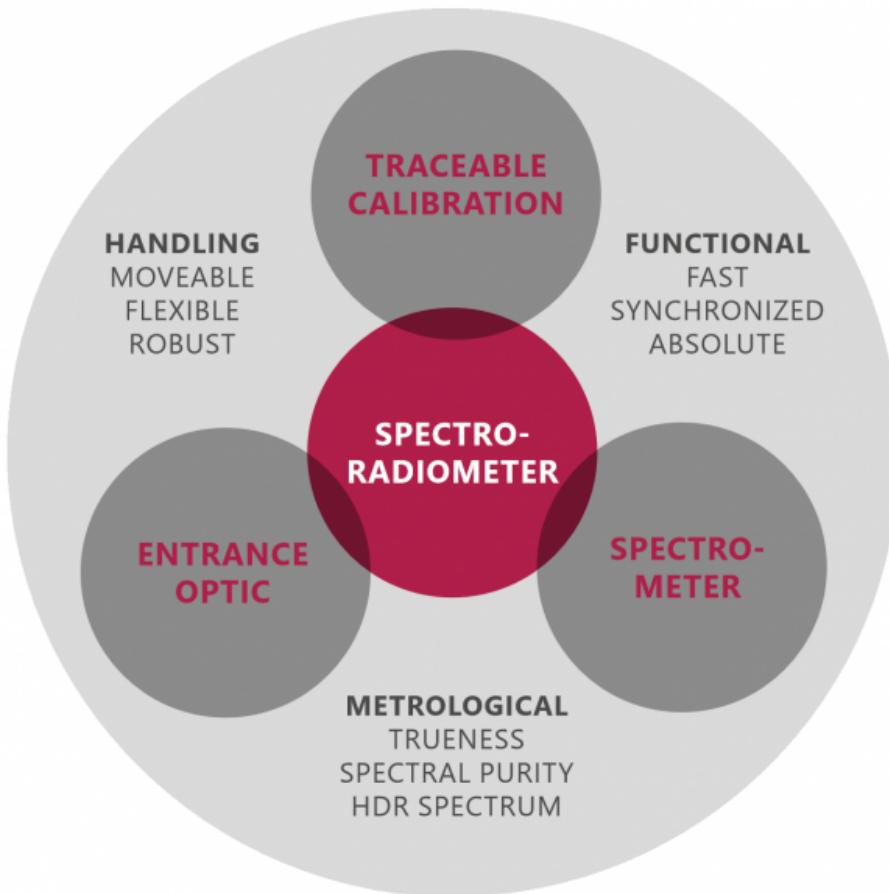

BTS2048-xx in Case




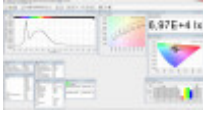
with SRT tube for 1°, 2°, 5°, 10° FOV

Downloads

| Type | Description | File-Type | Download |
|------------------|--|-----------|---|
| Brochure | Light measurement solutions for general and specialized lighting | pdf | https://www.gigahertz-optik.com/assets/Uploads/generallighting-broschuere-DINA4-hoch-v2.pdf |
| BTS2048 Brochure | Not Just Another Spectrometer | pdf | https://www.gigahertz-optik.com/assets/Uploads-v2/BTS2048-broschuere-DINA4-hoch-v2-WEB.pdf |

Configurable with

| Product Name | Product Image | Description | Go to product |
|--------------|---------------|--|---|
| CP-LG Series | | The CP-LG Series offers all kind of optical light guides. We have many different length, detector designs (inline, 90°, probes, etc.) available. Also customized designs are possible. | https://www.gigahertz-optik.com/en-us/product/cp-lg-series/ |

| Product Name | Product Image | Description | Go to product |
|---------------|---|--|---|
| S-SDK-BTS2048 |  | Software Development Kit for BTS2048 variants. | https://www.gigahertz-optik.com/en-us/product/s-sdk-bts2048/ |
| S-BTS2048 |  | Application software for BTS2048 variants. | https://www.gigahertz-optik.com/en-us/product/s-bts2048/ |

Contact, Calibration, Service & Support

We are known worldwide for excellent technical consulting and after sales support. Contact us to find together the best solution for you. Our services:

- Technical Consulting & Sales
- After-Sales Support
- Calibrations & Re-Calibrations ([ISO/IEC 17025 Calibration Services](#), [factory calibration](#), [Calibration of Third-Party Products](#))
- Repairs & Updates
- OEM & Feasibility Consulting of Customized Solutions

[Send us your inquiry](#) or contact us by phone or e-mail. We would welcome your feedback too or review us on [Google](#).

Gigahertz Optik GmbH (Headquarter)

Tel.: +49 (0)8193-93700-0
Fax: +49 (0)8193-93700-50
info@gigahertz-optik.de

An der Kaelberweide 12
82299 Tuerkenfeld, Germany

Gigahertz-Optik, Inc. (US office)

Phone: +1-978-462-1818
info-us@gigahertz-optik.com

Boston North Technology Park
Bldg B - Ste 205
Amesbury, MA 01913 USA